
Quick Start Guide


PRIMING UNIT (FIRST TIME USE ONLY)

Fill reservoir with warm water and replace it firmly on the base.

Insert tip firmly until it clicks into place.

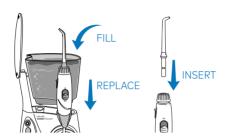


Turn pressure to high, point tip into sink, turn power ON and water switch to ON until water flows. Turn power and water OFF.

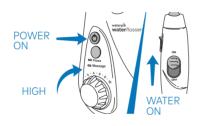

WATER FLOSS

With power OFF and water switch in OFF position, place tip in mouth. Turn power ON. Turn water switch to ON.

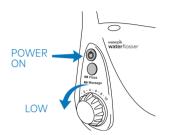

Lean over sink. Close lips enough to prevent splashing, while still allowing water to flow from mouth into sink



NEED HELP?

Visit **Waterpik.com/WelcomeWF**To watch a How To Video and find additional product support.

Call Us! 1-800-525-2774

Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: PLACE TIP
IN YOUR MOUTH BEFORE
TURNING WATER ON.

